

**BOWA**  
EINFACH SICHER

BOWA-electronic GmbH & Co. KG  
Heinrich-Hertz-Strasse 4 – 10  
72810 Gomaringen | Germany

Phone +49 (0) 7072-6002-0  
Fax +49 (0) 7072-6002-33  
info@bowa.de | bowa-medical.com




# ERGO 330R

## COAGULATION AND CUTTING IN A SINGLE INSTRUMENT



**BOWA**  
EINFACH SICHER

## The grasper with a unique added feature

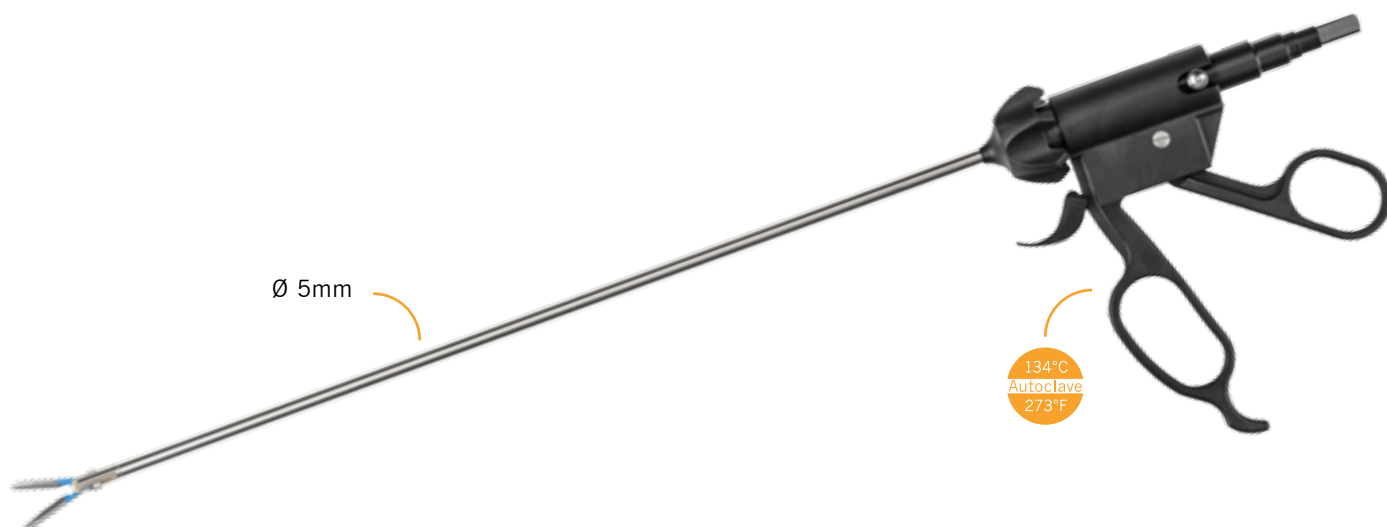

ERGO 330R is the intelligent combined bipolar instrument for laparoscopy from BOWA.

It allows tissue preparation, bipolar coagulation and mechanical tissue dissection without the complication of changing instruments. The ERGO 330R is a reusable instrument with a single use replaceable blade, for consistently good cutting performance with each operation. Different working lengths allow it to be used in various types of application.

The ERGO bipolar instrument is available with a standard flat plug connector and can be used with all common electro-surgical equipment.

In order to meet the constantly increasing requirements of modern hospitals for reprocessing and quality assurance, BOWA offers a COMFORT version with plug & play function and integrated service lifetime monitoring.



**The jaw insert** is serrated and enables secure gripping and holding of tissue. The tissue held by the device is cut centrally.





### **The BOWA COMFORT SYSTEM**

Immediate deployment by automatic parameter pre-selection on your ARC generator, for safe device settings and hence increased patient safety.

# Ordering Information



## Coagulation forceps with cutting blade Ø 5mm

| Grip variant   | Jaw insert          |                     | Shaft               |                     | Sliding rod         |                     | Blade holder        |                     | Complete Instrument consisting of:<br>- Handle<br>- Jaw insert<br>- Shaft<br>- Sliding rod<br>- Blade holder<br>- 1x743-999 Blade (qty. 20) |
|--|---------------------|---------------------|---------------------|---------------------|---------------------|---------------------|---------------------|---------------------|---|
|  | WL 340mm<br>743-334 | WL 250mm<br>743-325 | WL 340mm<br>743-234 | WL 250mm<br>743-225 | WL 340mm<br>743-434 | WL 250mm<br>743-425 | WL 340mm<br>743-534 | WL 250mm<br>743-525 |   |
|    |                     | ●                   |                     | ●                   |                     | ●                   |                     | ●                   | 743-853   |
| ERGO 330R handle,<br>COMFORT<br>4.5m cable<br><br>743-100                          | ●                   |                     | ●                   |                     | ●                   |                     | ●                   |                     | 743-851   |
|  |                     | ●                   |                     | ●                   |                     | ●                   |                     | ●                   | 743-854   |
| ERGO 330R handle<br><br>743-101  | ●                   |                     | ●                   |                     | ●                   |                     | ●                   |                     | 743-852   |

### Sterile single-use blade

For use with ERGO 330R coagulation forceps with cutting blade

|   |                                    |         |
|---|------------------------------------|---------|
|  | ERGO 330R, blade, sterile (qty 20) | 743-999 |
|---|------------------------------------|---------|