


MetraBag

The FDA (U.S. Food and Drug Administration) recommends the use of a tissue containment system for morcellation during hysterectomy and myomectomy.




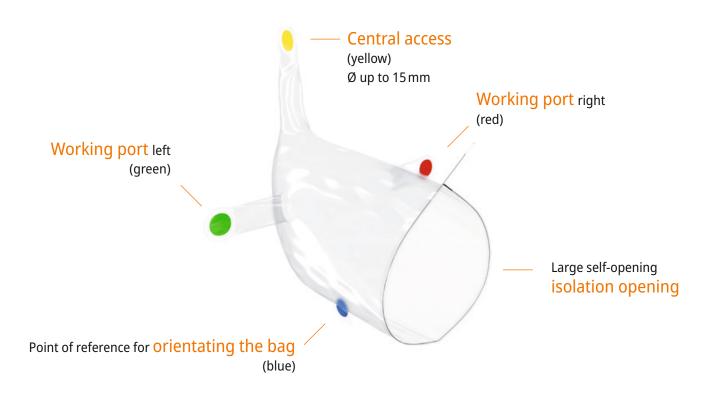

Intraoperative unrolling of the MetraBag

- MetraBag prevents residual tissue from being left behind in the abdomen without time-consuming suctioning.
- The variable working ports allow you to morcellate in a familiar way.
- The transparent MetraBag material provides a familiar working environment by maintaining visibility of surrounding structures.
- It is possible to introduce the MetraBag via a 10 mm trocar.



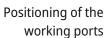



The MetraBag has been on the market for several years and has already been used thousands of times.

ERGO 300 Morcellator – precise, ergonomic, flexible

A powerful system for morcellation during hysterectomy and myomectomy.

BOWA MEDICAL offers the ERGO 300 with working diameters of 12 mm and 15 mm.

Vario foot switch with 3 m cable

Precise control of the cutting speed

Gear and sealing unit

Torque with uniform, quiet running



Safety

Deactivation of the blade for safe insertion into the operation area Handle for ergonomic control during the morcellation process

- Powerful morcellation system can be used for your surgical morcellation
- The ERGO 300 morcellation system is reusable and therefore economical.
- Complete system set up in the sterile area quick and easy.



Powerful micro motor with 3 m cable, autoclavable



The BOWA Morcellator can be used with the MetraBag tissue containment system.



Electrical loop for efficient removal of the cervix



Ordering information	REF
----------------------	-----

MetraBag		905-710
consisting of: Morcellation bag, insertion aid, 3 drawstrings sterile, 10 pieces		
Morcellator ERGO 300, set		905-000
Morcellator ERGO 300, set consisting of:		
Control unit ERGO 300 (TCM 3000 BL Morcellator)		905-001
Micro electric motor ERGO 300 with 3 m cable, autoclavable		905-002
Vario foot switch with 3 m cable		905-003
Gear unit with seals, complete		905-004
Handle		905-005
Obturator, Ø 15 mm		905-015
Protective tube, Ø 15 mm		905-215
Spray tube attachment		905-030
Spare part set, consisting of	Ø 12 mm	Ø 15 mm
Cutting tube		
Membrane seal Pack of 5	905-622	905-625
Roof-shaped seal		
Morcellator accessories	Ø 12 mm	Ø 15 mm
Obturator	905-012	905-015
Protective tube	905-212	905-215
Tenaculum forceps	905-412	
Claw grasper	905-415	
Morcellator seal with incision / 10 pieces	905-009	
Morcellator seal Ø 12–20 mm / 10 pieces	905-512	
Lubrication spray	905-032	
MetraLOOP		
MetraLOOP®, inner shaft, outer shaft		520-110
MetraLOOP®, COMFORT, inner shaft, outer shaft, cable 4.5 m	1	520-115
MetraLOOP® spare loops	Ø 100 mm	Ø 175 mm
MetraLOOP®, spare loops, single use, sterile / 10 pieces	520-117	520-118
MetraLOOP® cable		
BOWA, Erbe International, Martin International, Valleylab, Conmed (4.5 m)		370-050
Erbe T series and others with socket diameter of 4 mm (4.5 m)		280-035
Erbe VIO/ICC/ACC (4.5 m)		101-060



BOWA MEDICAL

BOWA-electronic GmbH & Co. KG

Heinrich-Hertz-Strasse 4–10 72810 Gomaringen, Germany

Phone +49 7072 6002-0

bowa-medical.com info@bowa-medical.com

BOWA MEDICAL provides the tissue morcellation solution for your gynaecological procedures

Contact your authorised BOWA medical device consultant now.

support@bowa-medical.com



EXPERTISE ensures SAFETY

Discover what the BOWA ACADEMY has to offer now




Checklist on HF surgery



Online course and certificate



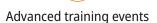



Live and on-demand webinars

bowa-medical.com/academy